Switch Adapted Water Gun




USER QUICK GUIDE

Introduction

This is a switch-adapted water gun that is compatible with 3.5mm switches. The intention of this device is to assist someone with limited hand function who might be interested in using a water gun.

Features

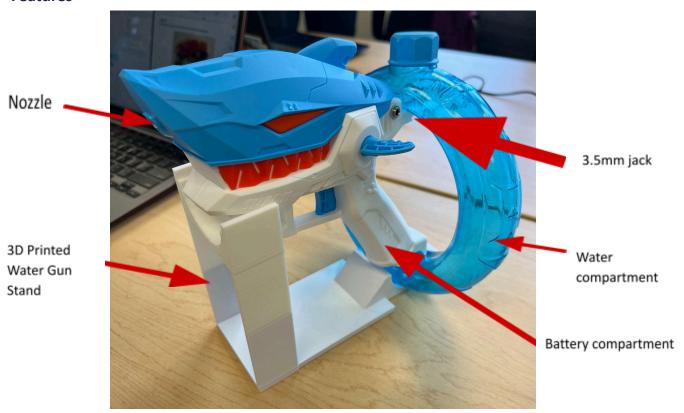

Usage

The Switch Adapted Water Gun can be used by filling up the water bottle with water and plugging in a 3.5mm switch into the side of the water gun. The water gun can be activated by either pressing the switch or using the trigger included with the water gun.

Compatibility

The Switch Adapted Water Gun should be compatible with any 3.5mm switch.

Switch Adapted Water Gun



USER QUICK GUIDE

Specifications

Length (mm)	Width (mm)	Height (mm)
177.8	50.8	177.8

Cleaning

The device can be cleaned using a damp cloth to wipe it down. The water bottle at the back of the water gun can be removed and emptied when not in use.

Changing the Batteries

The battery pack can be accessed on the side of the handle of the water gun by removing the screw to the compartment. Inside the compartment is a translucent container. Open the container to find the battery pack attached to the wires by a black plastic clip. Disconnect the battery pack by unclipping it and replace it with the new battery. The device uses a 3.7V rechargeable lithium battery.



